

REPORTED July 11, 2026

CLIENT

VERIFY Peptidetestcanada.ca/verify

EMAIL Info@peptidetestcanada.ca

ORDER # Matteo



END CUSTOMER

NAME Nexro Labs

EMAIL Nexrolab@gmail.com

PHONE —

ADDRESS —

§ 01 Sample Information

FOR RESEARCH PURPOSES ONLY

Sample ID: **BLD-025**

Compound: **GLOW**

Report No.: **COA-2026-SC-00214**

Lot Number: **COA-2026-SC-00214**

Total Mass: **151.5 mg**

Appearance: **Lyophilized powder**

Received: —

Analyzed: **07/11/2026**

Components: **3 peptides**

Total Claim: **70.00 mg**

Total Net: **81.69 mg**

SAMPLE PHOTOGRAPH




§ 02 Analytical Results

Qualitative and Quantitative chemical analysis by Liquid Chromatography Tandem Mass Spectrometry (LC-MS/MS)

PEPTIDE	LABEL CLAIM	NET CONTENT	CHROM. PURITY	LABEL ACC.	STATUS
GHK-Cu	50.00 mg	56.56 mg	99.29 %	113.1 %	CONFIRMED
BPC-157	10.00 mg	13.86 mg	99.44 %	138.6 %	CONFIRMED
TB-500	10.00 mg	11.27 mg	99.49 %	112.7 %	CONFIRMED
TOTAL NET CONTENT	70.00 mg	81.69 mg	—	116.7 %	

Chromatographic purity is reported as 100 % when no other peak is detected; if any additional peak is detected, the purity will vary from 100 %.

BLEND COMPOSITION share of total net peptide content

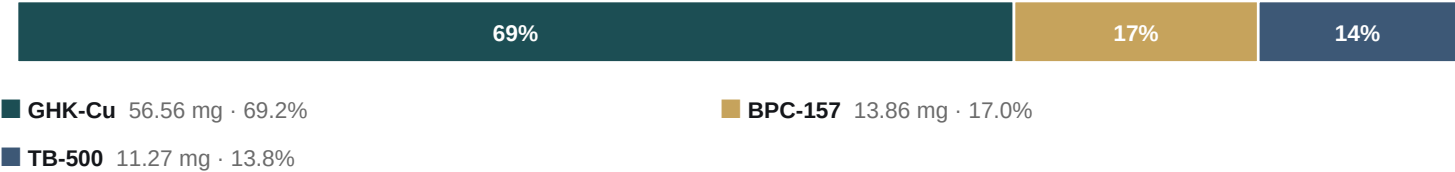


PRECURSOR [M+H]⁺

341.20 Da

Quantifier transition ● MRM

341.20 > 110.10 CE 15 V

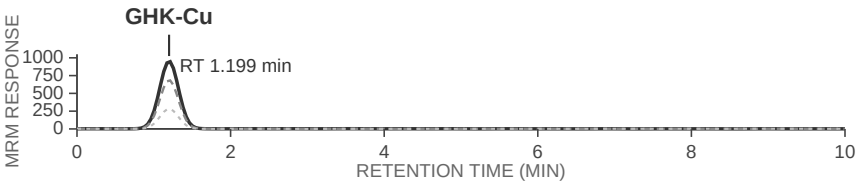

ROLE	TRANSITION	CE
● Quantifier	341.2 > 110.1	15 V
■ Qualifier	341.2 > 156.1	25 V


Agilent Poroshell 120 EC-C18 · 2.1×100 mm · 2.7 μm · Solvent A: 0.1 % formic acid in H₂O · Solvent B: 0.1 % formic acid in Acetonitrile

PRECURSOR [M+2H]²⁺

710.40 Da

Quantifier transition ● MRM

710.40 > 70.10 CE 22 V

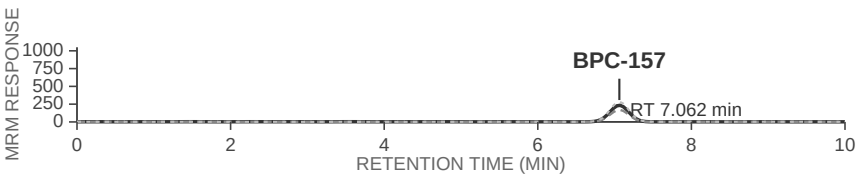

ROLE	TRANSITION	CE
● Quantifier	710.4 > 70.1	22 V
■ Qualifier	710.4 > 251.1	18 V


Agilent Poroshell 120 EC-C18 · 2.1×100 mm · 2.7 μm · Solvent A: 0.1 % formic acid in H₂O · Solvent B: 0.1 % formic acid in Acetonitrile

PRECURSOR [M+2H]²⁺

445.30 Da

Quantifier transition ● MRM

445.30 > 241.20 CE 12 V

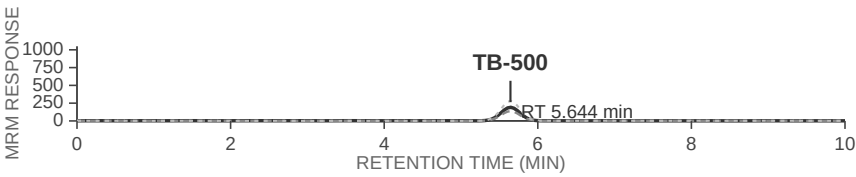

ROLE	TRANSITION	CE
● Quantifier	445.3 > 241.2	12 V
■ Qualifier	445.3 > 369.3	15 V

Agilent Poroshell 120 EC-C18 · 2.1×100 mm · 2.7 μm · Solvent A: 0.1 % formic acid in H₂O · Solvent B: 0.1 % formic acid in Acetonitrile

DISCLAIMER · RESEARCH USE ONLY

Results relate **only to the sample as received** and the tests performed — not to any other unit, lot, or batch. **For research use only:** not a certificate of quality, safety, efficacy, or regulatory compliance, and not for human or veterinary administration, clinical or diagnostic use, or label claims. SideChain Analytics does not manufacture, source, or sell the material and makes **no warranty, express or implied**; its liability shall not exceed the testing fee. This report may not be reproduced except in full.